

Technical Product Data Sheet											
<b>Tech. Sheet #</b>	003	<b>Version</b>	1.4	<b>Created:</b>	April 12 <sup>th</sup> 2016	<b>by</b>	R.Peralta	<b>Last updated</b>	October 6 <sup>th</sup> 2017	<b>by</b>	P.Cáceres
<b>Product Name</b>	Glutamate dehydrogenase (GDH)		<b>Product Code</b>	enz_gdh_003	<b>Current Dev. Phase</b>		Finished				
Core information											
<b>Product Type</b>	Lyophilized enzyme	<b>Producing microorganism</b>		<i>Escherichia coli</i> (recombinant)	<b>Microorg. code</b>	BL21	<b>Origin</b>	Thermophilic Bacteria			
<b>EC Number</b>	1.4.1.2			<b>CAS-No.</b>			9001-46-1				
<b>Product Description</b>	Oxidoreductase enzyme which relates carbon and nitrogen metabolism, catalyzing the reduction of $\alpha$ -ketoglutarate and ammonia to L- glutamate .										
<b>Temp Range °C</b>	20-70°C	<b>Opt. temp °C</b>	50°C	<b>Thermo stability</b>	Keeps more than 85% of its activity after 8 hours of exposure at 50°C	<b>pH range</b>	7.0-8.5	<b>Opt. pH</b>	8.0		
<b>Substrate</b>	$\alpha$ -ketoglutarate, NADH and ammonia										
<b>Products</b>	NAD <sup>+</sup> , glutamate, H <sub>2</sub> O										
<b>Reaction</b>	$\alpha$ -ketoglutarate + NADH + NH <sub>4</sub> <sup>+</sup> → Glutamate + NAD <sup>+</sup> + H <sub>2</sub> O										
<b>U (Unit definition)</b>	One unit is defined as the conversion of 1 $\mu$ mol of $\alpha$ -ketoglutarate into glutamate, in 1 minute at 50°C at pH 8.0.										
<b>Specific Activity</b>	≥ 90 U/mg protein										
<b>Protein concentration</b>	≥ 13% (w/w)										
<b>Molecular mass</b>	≈ 270 kDa			<b>Number of subunits</b>	Homohexameric (≈ 45 kDa subunit)						
<b>Substrate chirality</b>	No data available										
<b>Product chirality</b>	No data available										
<b>Alternative substrates</b>	No determined										
<b>Form</b>	Lyophilized powder										
<b>Other components</b>	0.05M Tris Buffer and 0.5 M NaCl (before lyophilizing)										
<b>Storage temperature</b>	-20°C										
<b>Stability</b>	At -20°C, it maintains the reported activity (≥ 90 U/mg) at least for 14 months										
<b>Shipping conditions</b>	Inside a styrofoambox with icepacks										

pH dependence

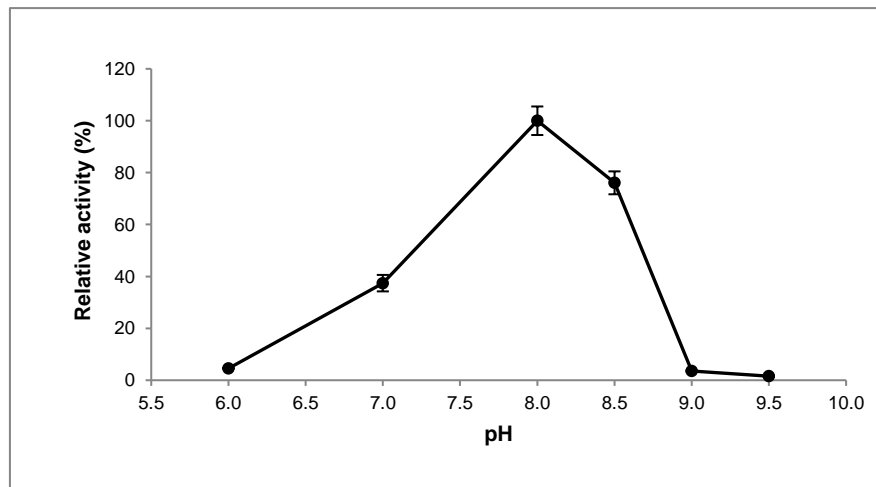


Fig 1. pH dependence of the rec GDH (enz\_gdh\_003). Activity was measured by monitoring pH from 6.0 to 9.5 at 50°C.

Temperature dependence

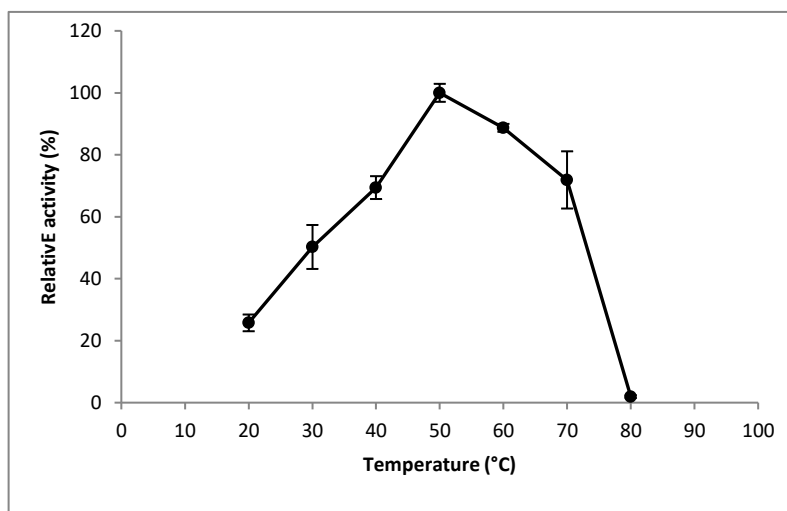
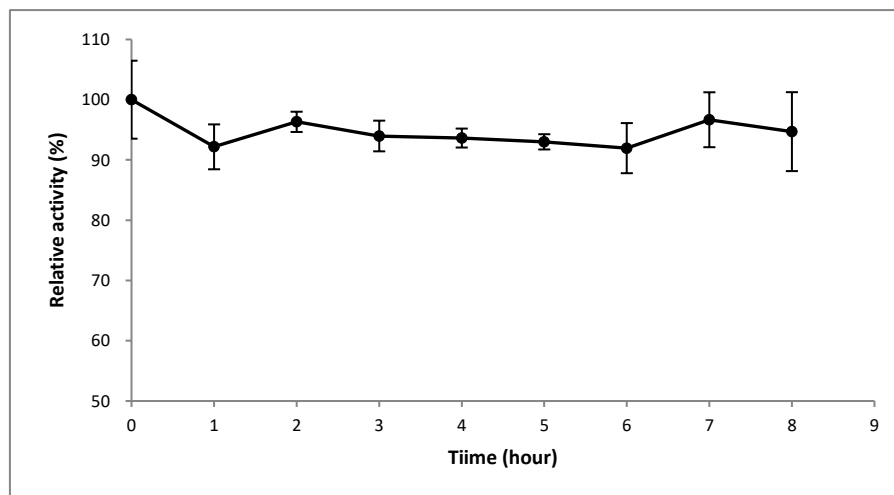


Fig 2. Temperature dependence of the rec GDH (enz\_gdh\_003). Activity was measured by monitoring temperature from 20 to 80° C in 50 mM Tris-HCl buffer (pH 8.0)

**Thermo-stability**



**Fig 3.** Thermostability of the rec GDH (enz\_gdh\_003). Activity was measured by monitoring at 50° C using 50 mM Tris-HCl buffer (pH 8.0)

**Scientific and Technical References**

1. D P. Hornby, M J. Aitchison, P C.Engel. (1984).The kinetic mechanism of ox liver glutamate dehydrogenase in the presence of the allosteric effector ADP. The oxidative deamination of L-glutamate. Biochemical Journal 1984-10-01
2. M. Amenábar, J. Blamey. (2011). Purification and characterization of a thermostable glutamate dehydrogenase from a thermophilic bacterium isolated from a sterilization drying oven. Biochemistry and Molecular Biology Reports. 2, 91-95.
3. J. DiRuggiero, F. Robb, R. Jagus, H. Klump, K. Borges, M. Kessel, X. Mai, M. Adams. (1993). Characterization, cloning, and in Vitro expression of the extremely thermostable glutamate dehydrogenase from the hyperthermophilic archaeon, ES4. Journal of Biological Chemistry. 268,17767-17774.

**Use, Security, Handling and Disclaimer**

<b>Use</b>	The product specified in this document has not been assessed for Human or Animal Consumption and is therefore <b>Not for Human or Animal Consumption</b> .
<b>Security</b>	The product specified in this document must be manipulated and used in a safe manner, in compliance with all applicable laws and regulations.
<b>Handling</b>	The product specified in this document must be manipulated by trained personnel, acquainted with laboratory security practices and regulations.
<b>Disclaimer:</b>	Any Product specified in this document is understood to be experimental in nature and may have hazardous properties. Swissaustral USA LLC MAKES NO REPRESENTATIONS AND EXTENDS NO WARRANTIES OF ANY KIND, EITHER EXPRESSED OR IMPLIED. THERE ARE NO EXPRESS OR IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OR THAT THE USE OF THE MATERIAL WILL NOT INFRINGE ANY PATENT, COPYRIGHT, TRADEMARK, OR OTHER PROPRIETARY RIGHTS. The user of the Product specified in this document assumes all liability for damages which may arise from its use, storage or disposal of this Product. Swissaustral USA LLC will not be liable to the user for any loss, claim or demand made by the user, or made against the user by any other party, due to or arising from the use of the Product by the user. Product specifications may change without notice.

Swissaustral USA LLC,  
111 Riverbend Rd #271, Athens, Georgia 30602, USA. Phone 706-206-8984  
e-mail: business@swissaustral.com